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4.9.1
 4.9.2

STEAMPlus™ Sterilizer Test Pack

For use in both gravity displacement and prevacuum steam sterilizers, the STEAMPlus™ Sterilizer Test Pack is equivalent in performance to an AAMI BI Test Pack (PCD).

Description

Designed to routinely challenge the steam sterilization process in healthcare facilities, each test pack contains a STEAMPlus™ Sterilization Integrator. Sterilizer test pack provides immediate read-out results and is "equivalent in performance to a AAMI BI Test Pack". When proper sterilization parameters have been met, the dark bar on the STEAMPlus Sterilization Integrator will enter the blue SAFE area.

Frequency of Use

For the highest level of sterility assurance, use a STEAMPlus Sterilizer Test Pack in every wrapped non-implant load.

Features & Benefits


- STEAMPlus™ Sterilization Integrator meets class 5 performance claims of ST-60 and ISO 11140-1:2005
- Immediate read-out results that requires no incubation
- Equivalent in performance to an AAMI BI Test Pack (PCD)
- Less expensive than a BI
- Monitors both gravity displacement and prevacuum sterilizers
- Load record card provides permanent, physical proof of sterilization conditions

Reference

ANSI/AAMI ST79. Comprehensive guide to steam sterilization and sterility assurance in health care facilities.

ORDERING INFORMATION

Cat. No.	Description	Unit of Measure
STP-025	STEAMPlus™ Steam Sterilizer Test Pack	25/cs

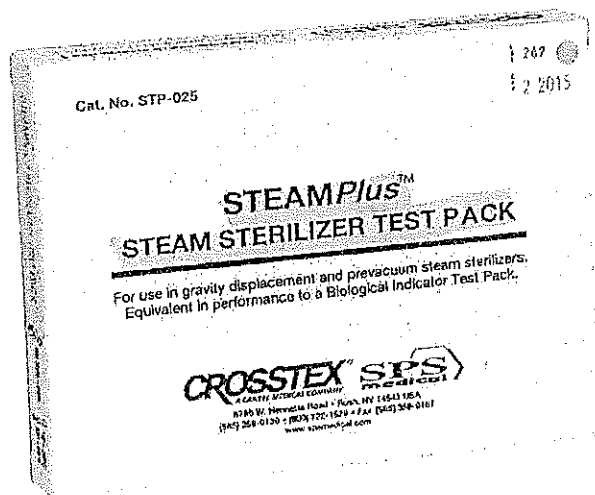STEAMPlus™ SPS medical
 STERILIZER TEST PACK LOAD RECORD CARD

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STEAMPlus™ STERILIZATION INTEGRATOR

BIOLOGICAL EQUIVALENT STERILIZATION TEST

If the dark bar has not reached the blue SAFE area, DO NOT release the load. Proper sterilization has NOT been achieved and the items in the load should be reprocessed. Prior to reprocessing, be sure to check the sterilizer for proper use and function. Maintain this load record card as physical proof of sterilization for this load.

Actual size is 4.5" x 6"

STEAMPlus is a trademark of SPSmedical Supply Corp.

SPSmedical Supply Corp. • 6789 W. Henrietta Road • Rush, NY 14543 USA
 (585) 359-0130 • (800) 722-1529 • Fax: (585) 359-0167 • E-mail: info@spsmedical.com
 Visit our website: www.spsmedical.com

Signature

STEAMPlus™ Type 5 Steam Sterilization Integrator

REF: SSI-100 / SSI-1000

INSTRUCTIONS FOR USE

Indication for Use:

The integrating indicator is designed to chemically react over time with the critical parameters of a steam sterilization cycle within a specified tolerance. The integrating indicator strip is intended to be placed in each pack, pouch, container, tray or other containment device to function as an independent monitor of critical parameters for the following sterilization cycles:

Standard Steam Sterilization Cycles		
Cycle Type	Exposure Temperature	Exposure Time
Dynamic Air Removal	270°F (132°C)	4 minutes
	270°F (132°C)	5 minutes
	270°F (132°C)	6 minutes
	270°F (132°C)	7 minutes
	270°F (132°C)	8 minutes
	270°F (132°C)	9 minutes
	270°F (132°C)	10 minutes
	273°F (134°C)	4 minutes
	275°F (135°C)	3 minutes
Gravity	250°F (121°C)	30 minutes
	270°F (132°C)	15 minutes
	275°F (135°C)	10 minutes
Immediate Use Steam Sterilization (IUSS) Cycles		
Dynamic Air Removal	270°F (132°C)	4 minutes
	275°F (135°C)	3 minutes
Gravity	270°F (132°C)	3 minutes
	270°F (132°C)	10 minutes
	275°F (135°C)	3 minutes
	275°F (135°C)	10 minutes

Device Description/Performance Characteristics:

The STEAMPlus™ Integrating Indicator is used in each pack to be sterilized to indicate exposure to specified critical parameters of steam sterilization. It meets the performance specifications for a Type 5 integrating indicator as defined in ANSI/AAMI/ISO 11140-1:2014.

Each strip consists of aluminum foil base adhered to label face stock that has been coated with a polypropylene film. A steam sensitive chemical containing a dark colored dye is held within

a well located at one end of the foil. A wicking strip leads from the well to the opposite side of the strip and during the steam sterilization cycle, steam penetrates the polypropylene film and comes in contact with the steam sensitive chemical. The steam sensitive chemical requires all the critical parameters of steam (Steam, Temperature and Time) to melt and wick along the filter paper carrying the dark dye with it. This is made visible through the window. The integrator strip indicates that the specified critical parameters of steam sterilization have been achieved when the dark colored steam sensitive chemical containing dye reaches the "PASS" window.

The STEAMPlus™ Integrating Indicator and its packaging are manufactured in the absence of lead or other heavy metals. It is also manufactured without natural rubber latex and dry natural rubber.

Instructions for Use:

Wrapped items, Rigid Sterilization Containers and Sterilization Pouches

1. Place an integrator strip(s) inside each pack, pouch, container, tray or other containment device, in the area(s) of greatest challenge to the sterilant but not directly beneath an item of great metal mass.

NOTE: Place an indicator strip on each level of multi-level trays or cassettes.

NOTE: If the packaging is opaque, use an external process indicator.

2. Follow the containment devices instructions for use to prepare for steam sterilization.

3. Upon completion of the cycle, follow department procedure for load release.

NOTE: Confirmation of the integrator strip result should take place at time of pack use.

4. Prior to use of the items within the pack, aseptically remove the integrator strip(s) and observe the location of the dark bar within the integrators window.

a.) The integrator demonstrates passing results when the dark bar has completely travelled through the "FAIL" area or window and has entered the "PASS" area.

NOTE: Variability in the color of the dark bar is normal and acceptable as long as the dark color is visible in the "PASS" area.

b.) The integrator demonstrates failing results when the dark bar is not visible in the "PASS" area. Follow departmental procedures for investigating suspected sterilization failures.

5. See Reference Guide for examples of pass and fail results.

SteamPlus™ Type 5 Steam Sterilization Integrator


Introduction

The SteamPlus™ Sterilization Integrator is an advanced technology chemical indicator that provides a simple, accurate method of assuring that proper conditions for sterilization have been met during a sterilizer cycle. For use with all steam processes (gravity, prevacuum and IUSS), the SteamPlus Integrator has documented performance equivalent to a biological indicator. Strip lot number and expiration date are printed on front side for easy reading.

Description

The SteamPlus Sterilization Integrator is engineered to integrate the three (3) critical variables of sterilization: time, temperature, and saturated steam. This information is displayed in a precise, easy-to-read format. When the dark bar enters the blue SAFE area, sterilization criteria has been met. If the dark bar does not reach the blue SAFE area, proper conditions for sterilization have not been achieved and the processed pack or load should not be released. The SteamPlus Sterilization Integrator has been tested to the AAMI ST-60 and ISO 11140-1:2014 standards performance requirements; Note: The FDA does not recognize Type 3, 4 or 5 indicators as defined in ISO 11140-1:2014.


SteamPlus™ Type 5 Integrators		
Ref. #	Description	Quantity
SSI-100	4 in x 0.75 in Strip (10.2 cm x 1.9 cm)	100 Count
SSI-1000	4 in x 0.75 in Strip (10.2 cm x 1.9 cm)	1000 Count

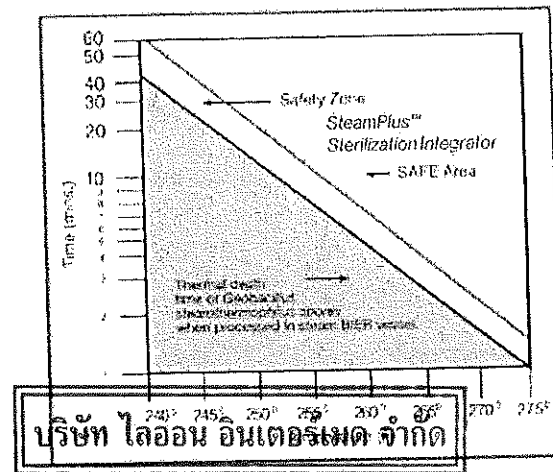How to Use the SteamPlus™ Steam Sterilization Integrator

1. Place a SteamPlus Integrator inside each pack or minimally in every load (i.e. middle of the load or inside the densest pack) to be sterilized.
2. Process the load according to the sterilizer manufacturer's instructions.
3. Release the load or pack if the dark bar on the SteamPlus Integrator has entered the blue SAFE area. If the dark bar has not entered into the blue SAFE area, DO NOT release the load or pack.

While the SteamPlus Integrator is FDA cleared as equivalent in performance to a biological indicator, it does not replace the need to run a weekly or preferably daily BI test. However, it can be used as part of the immediate load release criteria for non-implant loads. Implant loads should be quarantined until BI results are known.

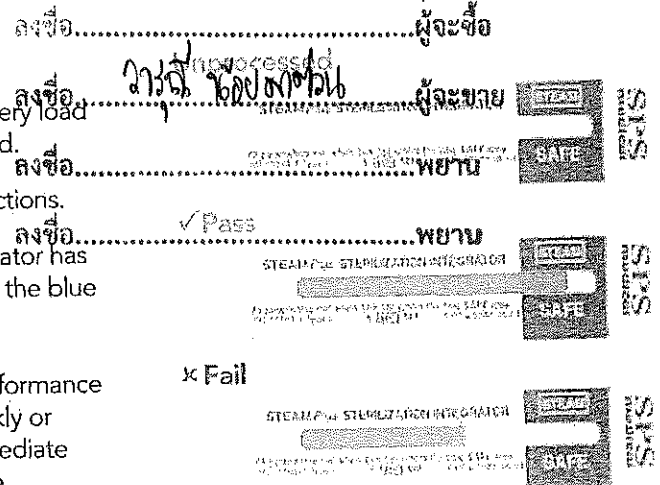All company and product names are trademarks or registered trademarks of their respective owners. All rights reserved. All trademarks are the property of their respective owners.

© 2021 Crosstex Medical Corporation, Inc. All rights reserved. All trademarks are the property of their respective owners.

STEAMPlus™ Class 5 Steam Sterilization Integrator

Introduction

The STEAMPlus Sterilization Integrator is an "advanced technology" medical device that provides a simple, accurate method of assuring that proper conditions for sterilization have been met during a sterilizer cycle. Certified for use with all Steam processes (gravity, prevacuum and flash), the STEAMPlus has documented performance equivalent to a biological indicator.

Description

The STEAMPlus Sterilization Integrator is engineered to integrate the three (3) critical variables of sterilization: time, temperature, and saturated steam. This information is displayed in a precise, easy-to-read format. When the dark bar enters the blue SAFE area, sterilization criteria has been met. If the dark bar does not reach the blue SAFE area, proper conditions for sterilization have not been achieved and the processed items should not be released.

Some common causes of sterilization failure are:

- overloading the chamber
- air pockets within the sterilizer
- packs wrapped too tightly
- malfunction of the timing mechanism
- error in the temperature setting
- load configuration
- poor steam quality

Technical Design

The STEAMPlus Integrator is made of aluminum foil, which holds the temperature and steam-sensitive chemical which is designed to melt when subjected to a steam environment. The foil also acts as a moisture barrier against steam during the sterilization cycle.

As moisture penetrates the integrator, it lowers the melting point of the chemical. When melting occurs, the liquid chemical is soaked up by the paper wick and, as time elapses, moves along the scale. The rate of melt produced is a function of both the moisture-vapor transmission rate of the cover film and the melting point depression of the chemical. The combination of these two factors provides a rate of melt at various temperatures which parallels the thermal death time of *Geobacillus stearothermophilus* spores as illustrated in Figure 1.

Note: Figure 1 represents sterilization cycle time only and does not include come-up or down time.

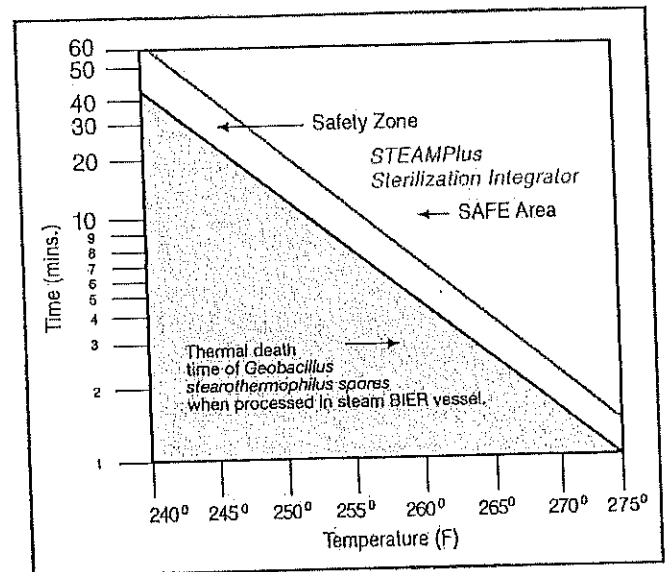
Figure 1

Performance Characteristics

The STEAMPlus Integrator was tested at various time and temperature intervals in saturated Steam in order to estimate the time required at each temperature for the STEAMPlus Integrator to produce a SAFE reading. Five temperature points between 245°F and 275°F were tested for at least five separate time periods. This testing resulted in the following mean (average) times to produce a SAFE reading for each of the five temperatures tested:

245°F/121°C	@ 21.0 minutes
250°F/121°C	@ 21.0 minutes
260°F/127°C	@ 7.1 minutes
270°F/132°C	@ 2.3 minutes
275°F/135°C	@ 1.5 minutes

The performance of the STEAMPlus Integrator in comparison with the performance of *Geobacillus stearothermophilus* spores, plus a slight safety margin, is shown above in Figure 1. This margin of safety, which spans the entire spectrum of normal Steam sterilization temperatures, provides added assurance that once the dark bar enters the blue SAFE area, sterilization has indeed occurred. The STEAMPlus Sterilization Integrator meets Class 5 performance claims of AAMI ST-60 and ISO 11140-1:2005. Note: The FDA does not recognize Class 3, 4 or 5 in ISO 11140-1:2005.

Quality Control

When tested in a Steam BIER (Biological Indicator-Evaluator Resistometer) vessel for 1 minute at 270°F (132°C), the STEAMPlus Integrator is designed to not register in the blue SAFE area. Conversely, when tested at 270°F (132°C) for 3.0 minutes, the STEAMPlus Integrator is designed to register in the blue SAFE area, thus signifying sterilization has been achieved. A production lot number & expiration date is printed on each device for traceability.

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Immediate Use Steam Sterilization:

1. Place the integrator strip in the center of each configuration (tray, protective case, flash approved rigid sterilization container) to undergo immediate use steam sterilization.
2. Follow the containment device's instructions for immediate use steam sterilization.
3. Upon completion of the sterilization cycle, review the printout for appropriate sterilization parameters. Remove the immediate use sterilized item(s).

CAUTION: THE CONTAINMENT DEVICE WILL BE HOT. USE THERMAL PROTECTIVE GLOVES.

4. Prior to use of the items within the pack, aseptically remove the integrator strip(s) and observe the location of the dark bar within the integrators window.

a.) The integrator demonstrates passing results when the dark bar has completely travelled through the "FAIL" area or window and has entered the "PASS" area.

b.) The integrator demonstrates failing results when the dark bar is not visible in the "PASS" area. Follow departmental procedures for investigating suspected sterilization failures.

5. See Reference Guide for examples of pass and fail results.

Storage Conditions:

Before use, store all STEAMPlus™ Integrating Indicators in their original closed zipper bag at room temperature. Do not store integrator strips near heat, moisture, oxidizing agents, acids/alkalis or cleaning/disinfecting agents.

The endpoint color of exposed integrator strips is stable for 36 months when stored away from light at room temperature.

Expiration Date:

Do not use beyond the expiration date printed on the integrator strip. The expiration date is printed starting with the four digit year followed by the two digit month ending with the two digit day.

Example: 20200509 expires the year 2020 on May 9th.

Performance Limitations:

- This product is for single use only. Do not reuse partially processed integrator strips.
- The STEAMPlus™ Integrating Indicator is part of a quality control system for sterility assurance. They cannot be used as a sole means for validating the sterilization process.
- All personnel involved must be trained in the proper use of the integrator strips and a system for reporting failures must be established.
- If there are any problems interpreting the results from the product, contact HuFriedyGroup or your local HuFriedyGroup representative.

Safety Precautions:


This device contains a steam sensitive chemical and indicator dye that is safe under normal conditions of use. This device is not intended for human consumption.

Disposal:


This product is manufactured in the absence of lead and other heavy metals and may be disposed of as regular waste.

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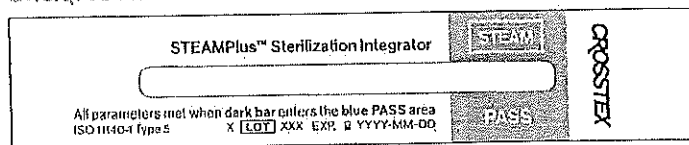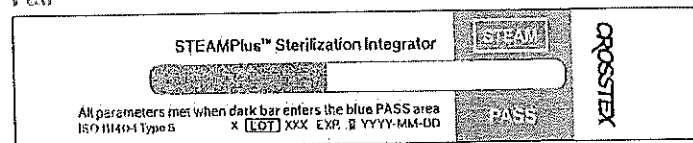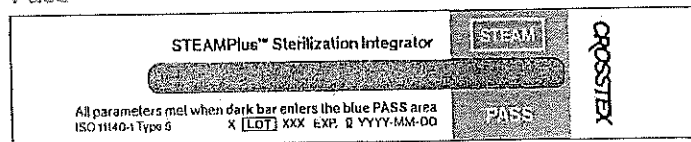
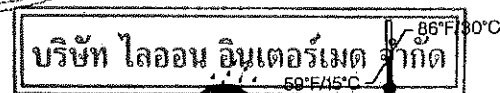
CROSSTEX

Manufactured for: Crosslex International, Inc.
6789 W. Henrietta Road, Rush, NY 14543 USA
585.359.0130 | 631.582.6777

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CE STERIS GmbH c/o BDO AG
Längelweg 116A
CH-8504 Biel-Bienne
SWITZERLAND

CE STERIS Solutions Limited
Chancery House, Rayns Way
Watermead Business Park, Syston
Leicester
LE17 1PF
United Kingdom

Made in USA
103501-HFU / 01/2023

In case of any serious incident that has occurred and involves our devices within the European Union, please report it to HuFriedy Mfg. Co., LLC and your national competent authority of the member state.
All company and product names are trademarks of HuFriedy Mfg. Co., LLC, its affiliates or related companies, unless otherwise noted.
©2023 HuFriedy Mfg. Co., LLC. All rights reserved.

HuFriedyGroup
The Best In Practice

4.9.2

4.9.2

REF STP-025

STEAMPlus™ STEAM STERILIZER TEST PACK

A Type 5 Integrator for use in gravity displacement and prevacuum steam sterilizers. Equivalent in performance to a Biological Indicator Test Pack.

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Frequency of Use

For the highest level of sterility assurance, use a STEAMPlus™ Sterilizer Test Pack in every non-implant load.

Instructions for Use

1. Place a Test Pack in the sterilizer on the lowest shelf over the drain along with a full load.
2. Process the load according to the sterilizer manufacturer's instructions.
3. After processing, remove the Test Pack and allow to cool before opening.
4. Confirm that the indicator on the Test Pack has turned brown and then remove the STEAMPlus™ Load Record Card from the center of the Pack.
5. If the dark bar has not entered into the blue SAFE area, the load must be reprocessed.
6. Complete the Load Record Card and retain as physical proof of sterilization.

Storage and Shelf Life

Store in controlled room temperature away from sterants. The STEAMPlus™ Sterilizer Test Pack has a two (2) year shelf life from date of manufacture.

All company and product names are trademarks of Crosstex International, Inc.
All other trademarks are the property of their respective owners.
©2001 Crosstex International, Inc. All rights reserved.

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
[Handwritten signature]

ProEZ AW Quad™

Concentrated Quadruple Enzyme Automatic Washer Detergent

Description:

ProEZ AW Quad is a four enzyme, non-foaming, neutral pH detergent containing powerful surfactants for use in automatic washers. Protease, amylase, lipase and cellulase enzymes work together with carefully selected surfactants to clean surgical instruments and endoscopes.

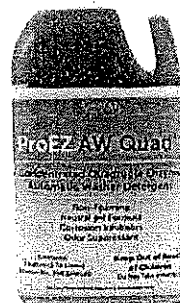- Four enzymes attack the entire range of surgery soils: blood, fats and starches.
- Verified enzyme activity at expanded temperature range from 60°-150°F (16°-65°C).
- Formula delivers the concentration of active agents needed for fast, consistent results. Starts working quickly - in as little as one minute in warm water.

Intended Use:

ProEZ AW Quad is specially formulated for use in high impingement automated washers and washer/sterilizers. ProEZ AW Quad is intended for cleaning medical instruments by dissolving and lifting organic soils in as little as 1 minute. It is specifically designed for cleaning instruments with difficult to access areas such as endoscope channels, box locks and serrated edges.


Product Specifications

Dilution Ratio	1/2 oz. to 1 oz. per gallon
pH	Neutral
Temperature of Use	ProEZ AW Quad has an expanded temperature range. It is effective at room temperature about 68°F (20°C). Use with warmer water at 90-110°F (32-43°C) will optimize the performance of the enzymes. This product may be used in automated cleaning equipment (ultrasonics, AEDs, instrument washers) from 90-150°F (32-65°C) according to equipment manufacturer's instructions
Odor	Fresh
Color	Blue to green
Shelf Life	2 years from date of manufacture

Effective in multiple applications:

- Presoak and manual sink cleaning
- Ultrasonic instrument cleaning devices
- Manual endoscope cleaning and automated cleaning cycles of AERs
- Clean and deodorize surgical suction units and lines
- Automated washers for main cycle

Premium enzymatic detergent provides the chemistry needed for difficult soils, complex instrumentation and high impingement automated cleaning equipment.

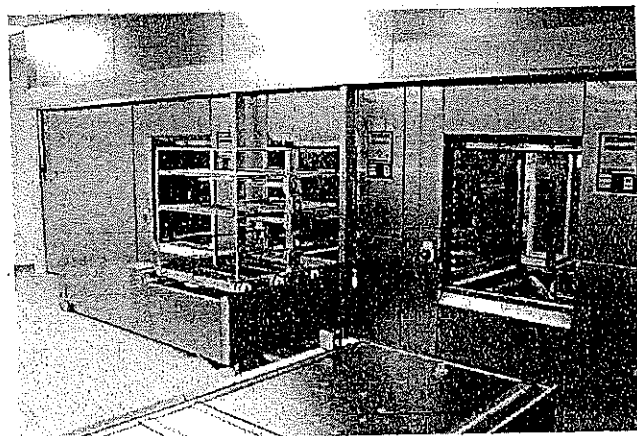Description	Case Qty	Order#
1 Gallon	4	PREZAW128-1
1 Gallon Re-Order w/o Pump	4	PREZAWN128-1
5 Gallon Drum	1	PREZAW050-1
15 Gallon Drum	1	PREZAW150-1
30 Gallon Drum	1	PREZAW300-1
55 Gallon Drum	1	PREZAW550-1

[Handwritten signature]

ProEZ AW Quad™ Concentrated Quadruple Enzyme Automatic Washer Detergent

WHAT IS PROEZ AW QUAD?

ProEZ AW Quad is a premium non-foaming neutral pH, multi-enzymatic detergent. Cleaning action is enhanced by a full range of enzymes needed to dissolve complex soils found on surgical instruments: protease acts on proteins, amylase acts on carbohydrates, lipase acts on fats, and cellulase acts on the indigestible fiber and biofilm in human soils. Enzyme action is enhanced by penetrating agents, chelating agents to improve cleaning action in hard water and a system of anti-corrosive agents.

ProEZ AW Quad is indicated for the most challenging gastrointestinal and orthopedic soils.

What are the recommendations for using ProEZ AW Quad?

- Use in the first or "enzyme" cycle of automated high impingement instrument washers.
- Use in ultrasonic instrument cleaning equipment.
- Use as pre-soak or holding solution to keep soils moist on contaminated surgical instruments and related equipment until further processing can take place.
- Use as a "flush and soak" solution for manual cleaning of all types of immersible endoscopes.
- Use in the cleaning cycle of automated endoscope reprocessors. The no-foam formula protects the pumps and flushing mechanisms in automatic reprocessors. Always refer to the manufacturer's recommendations.
- Use for manual cleaning of all types of immersible surgical instruments and equipment. The no-foam formula offers increased safety by allowing technicians to visualize instruments in the cleaning pan or sink.

How is ProEZ AW Quad different from ProEZ 1™ & ProEZ 2™?

- It is specifically formulated to be non-foaming and suitable for use in high impingement automated instrument washers.
- The premium formula delivers higher concentrations of enzymes supported by a chemical team of penetrating and rinsing agents. Because the components enhance and intensify enzyme cleaning action, the formula is suitable for the most challenging soils and instrumentation.

- It is the ideal solution for standardization because it can be used anywhere in any system.
- The non-foaming formula contributes to safety by allowing instruments to be visualized during manual processing.

What is the ideal dilution when using ProEZ AW Quad?

It is effective starting at a dilution ratio as low as 1/4 ounce per gallon of water. Use the recommendations of automated instrument washer manufacturers as a starting point. Washer sterilizers and European manufactured washers with low water pressure spray arms may require higher percentages of detergent per cycle than high impingement washers with strong spray action. Water quality is also a key factor in determining the appropriate dilution. "Hard" water inhibits cleaning by consuming some of the chemical action of detergents.

Water is rated as "hard" (i.e. higher concentrations of calcium carbonate and other minerals) at 251 ppm or higher. Water rated as "hard" will require higher concentrations of detergents to yield effective results. Water quality information is usually available from the local water department.

What is the recommended temperature range when using ProEZ AW Quad?

ProEZ AW Quad is a uniquely formulated multi-enzymatic detergent with expanded temperature capabilities.

- Manual Use: (pre-soak pans, flushing cannulated items, manual sink scrub) Dilute ProEZ AW Quad with water at a minimum of 68-110°F (20-43°C). In warm water ProEZ AW Quad works in as little as 1 minute making it the fastest acting enzymatic detergent on the market.
- Automated cleaning equipment: ProEZ AW Quad may be used in ultrasonic instrument cleaning equipment and automated endoscopic reprocessors (AERs) set from 90-120°F (32-49°C). The diluted cleaning solution should be changed more frequently when used at higher temperatures.
- Automated washers: ProEZ AW Quad is effective for the enzyme or pre-wash cycle at 90-120°F (32-49°C). It is tested and validated to provide effective cleaning in the main wash cycle up to 150°F (65°C).

Four Powerful Cleaning Enzymes:

- Protease, amylase, lipase and cellulase enzymes work together to quickly dissolve protein, starch and fatty soils in as little as one minute.
- Ideal formulation for your toughest cleaning challenges such as orthopedic and gastrointestinal soils.
- Thorough cleaning for complex instrument surfaces including cannulated areas, box locks and serrations.

Superior Cleaning System:

- Carefully selected surfactants reduce surface tension, improve penetration of cleaning agents and promote free rinsing of dissolved soils.

Multi-Use Detergent:

- Simplify purchasing and inventory. Standardize on one highly effective formula.
- Ideal for high impingement automated washer/disinfectors and washer/sterilizers.
- No foam formula also indicated for ultrasonic cleaning systems including automated washers for cannulated items.
- Fast, effective cleaning of tough gastrointestinal soils in automated endoscope washers.
- No foam formula improves visibility for increased safety during manual/pan soaking and flushing.

Cost Effective:

- ProEZ AW Quad is highly concentrated with effective dilution starting as low as 1/2 ounce per gallon. Adjust dilution based on soil load and water quality.

Product Capabilities:

ProEZ AW Quad rapidly digests protein, starch and fats, the principle components of surgical soils, and is non-corrosive and safe to use on all instruments.

It is a convenient 4-in-1 product which may be used:

1. As an instrument holding solution.
2. To clean instruments in manual systems.
3. In ultrasonic cleaners and scope reprocessors.
4. In automatic washers.

The non-foaming solution aids in the visualization of instruments when used in manual systems.

Use Precautions:


As with other enzymatic detergents, ProEZ AW Quad is an eye irritant and mild skin irritant. Wear eye protection and water resistant utility gloves when using the product.

Availability:

Stocked by MicroCare Medical distributors. Contact MicroCare Medical for a distributor near you or visit MicroCare.com/Certol.

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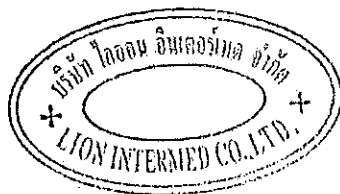[Handwritten signature]

MicroCare Medical
6120 East 58th Avenue
Commerce City, CO 80022
(800) 843-3343

Visit Us:
MicroCare.com/Certol

© 2022 MicroCare, LLC. All Rights Reserved. "MicroCare Medical", "ProEZ AW Quad" are trademarks or registered trademarks of MicroCare, LLC. The information set forth herein is based on data believed to be reliable. MicroCare Medical makes no warranties express or implied as to its accuracy and assumes no liability arising out of its use by others. This publication is not to be taken as a license to operate under nor infringe upon any patents not herein expressly described.

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MicroCare is a versatile non-foaming, neutral pH, four enzyme detergent for cleaning dental instruments and devices. It dissolves soils making it easier to clean complex areas where biofilm can form. Standardize on one product for multiple cleaning tasks: manual cleaning, ultrasonic instrument cleaning, manual and automated cleaning of endoscopes, automated instrument washer enzyme and manual cleaning of surgical or dental evacuation systems.

MicroCare's enzymatic (protease, amylase, lipase and cellulase) combine with a powerful oxidizing agent to break down and lift complex soils especially in lumens or other hard to reach areas. It is safe for use on orthopedic implants. The neutral pH formula with anti-corrosives is safe for use on instruments including aluminum, rubber, plastic, tubing and optics.

For more information call toll free 1-800-535-5053. Have the product container or label with you when you visit your Central Center or doctor or going for treatment.

MicroCare is committed to providing high quality products. For a copy of any MicroCare product literature call 1-800-843-3343.


MicroCare's ingredients: Propylene Glycol, Subtilisin, Alpha-Amylase, Borax Decahydrate, Sodium Lauryl Sulfate, Phosphoric Acid, Diethyl phthalate, Dioctylamine, Ethanol.

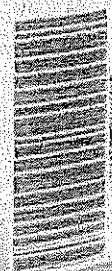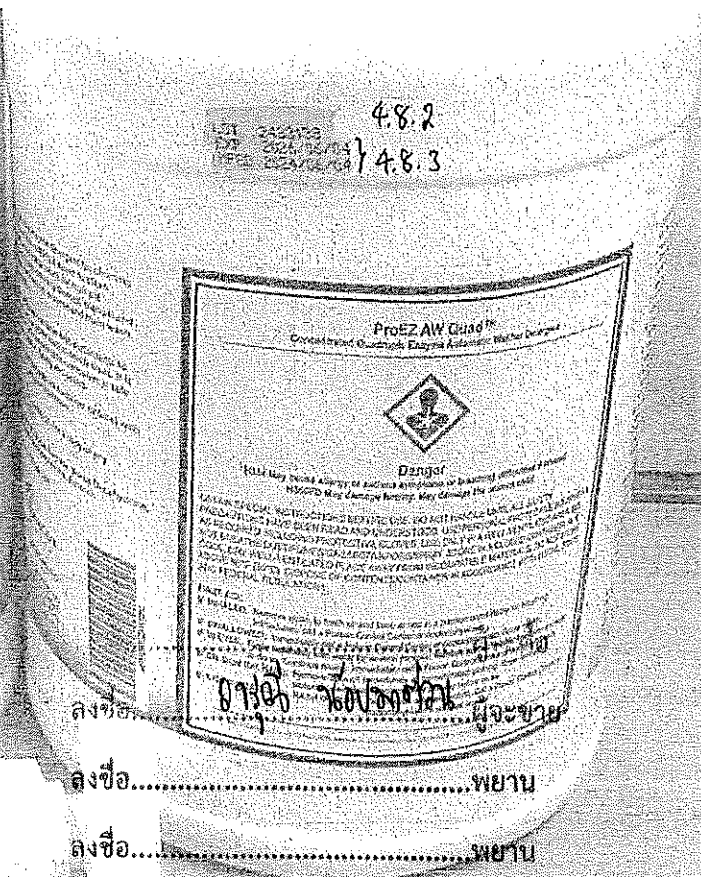
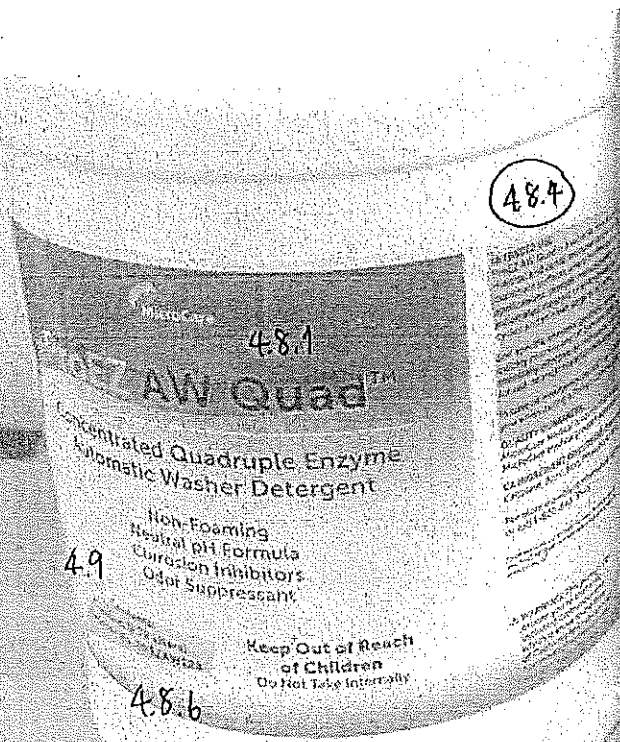
For more information visit www.MicroCare.com/eb-258

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